

ZPure™ Glass H₂O/HC filters

Compression fittings

Introduction

Pure gas is critical requirement in gas chromatography, spectroscopy, optics, lithography, and numerous other applications in manufacturing and analytical laboratories. The ZPure™ H₂O/HC filters remove H₂O and hydrocarbons (C5 and heavier) from dry air or inert gases.



Important



Warnings must be read carefully and understood. Improper use of this product can cause harm or death to personnel and damage to property!

- Wear eye protection and use caution when working with pressurized systems.
- Maximum pressure is 11 bar / 160 psi.
- When using a high-pressure gas source over-pressure protection must be provided.
- Do not open filter, even after use.
- Do not use with supply gas H₂O concentration greater than 0.1%. Contains molecular sieve which may become very hot in contact with liquid water. Do not use with liquids.
- Do not expose the protective outer plastic tube to solvents.
- Do not use without the protective outer plastic tube in place.

- Special precautions might be required when using hydrogen or hazardous gases. Consult local regulations and your company's safety procedures.
- Color Disclaimer: This product is intended to be yellow. Slight variations in shade may occur due to the manufacturing process and differences in individual color perception.

Instructions

Installation for compression fittings

1. Check the package contents for damage. Contact your supplier if any items are damaged or missing.
2. The ZPure filters are delivered in the active state, filled with helium. It is ready for use, but it is recommended that the trap be purged with the gas to be purified.
3. For the best connections, start with freshly-cut tubing in the gas line. 1/8" nuts should be turned 3/4 turn past finger tight. 1/4" nuts should be turned 1-1/4 turns past finger tight.
4. To connect in-line, first purge the gas line with supply gas. While maintaining a low purge rate of 10 to 30 cc/min, remove the end plugs from both ends of the trap, then attach the inlet end of the trap to the gas line.
5. While the end plugs are removed from the trap during installation, a small amount of air will diffuse into the adsorbent. As long as the trap is attached to the gas line within a few minutes, the loss of adsorbing capacity will be negligible.
6. Attach the outlet end of the trap.
7. While maintaining a low purge rate of 10 to 30 cc/min, purge the entire system until a total of 1 to 2 liters of gas has been flushed through.
8. Leak check all connections.

Specifications

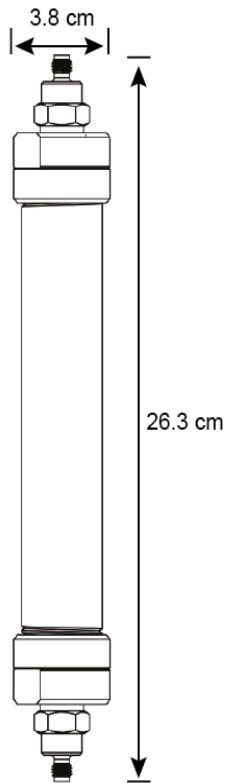


Figure 1. 46 cc

References



Recycling

Please contact Trajan for information regarding recycling purifier products.

Information and support

Visit www.trajanscimed.com or contact techsupport@trajanscimed.com

Specifications are subject to change without notice.

US AGENT	Trajan Scientific Americas Inc 2601 Technology Drive, Louisville, KY 40299, USA Tel: +1 502 491 6300
---------------------	--

Product Numbers

46 cc	
1/8" Brass	1/4" Brass
202256-B	202257-B
1/8" Stainless steel	1/4" Stainless steel
202256-SS	202257-SS